



TEL:4006-871-227 Web:www.ybio.net Email:shybio@126.com

YBA777Hu01 50 μ g
Recombinant Inducible T-Cell Co Stimulator (ICOS)
Organism Species: Homo sapiens (Human)

Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

10th Edition (Revised in Jan, 2014)

[PROPERTIES]

Residues: Phe31~Lys190

Tags: Two N-terminal Tags, His-tag and GST-tag

Accession: Q9Y6W8

Host: *E. coli*

Subcellular Location: Cell membrane; Single-pass type I membrane protein. Secreted.

Purity: >95%

Endotoxin Level: <1.0EU per 1 μ g (determined by the LAL method).

Formulation: Supplied as lyophilized form in 20mM Tris, 150mM NaCl, pH8.0, containing 1mM EDTA, 1mM DTT, 0.01% sarcosyl, 5% trehalose, and preservative.

Predicted isoelectric point: 7.6


Predicted Molecular Mass: 50.1kDa

Applications: SDS-PAGE; WB; ELISA; IP.

(May be suitable for use in other assays to be determined by the end user.)

[USAGE]

Reconstitute in ddH₂O.

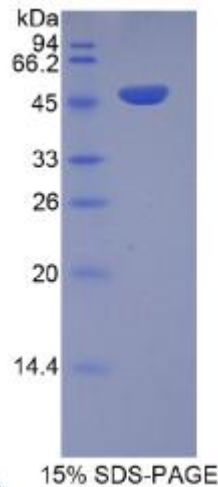



TEL:4006-871-227 Web:www.ybio.net Email:shybio@126.com

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate of the target protein. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. (Referring from China Biological Products Standard, which was calculated by the Arrhenius equation.) The loss of this protein is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCES]

The sequence of the target protein is listed below.

**FIFHNGGVQI LCKYPDIVQQ FKMQLLKGGQ ILCDLTKTKG SGNTVSIKSL KFCHSQLSNN
SVSFFLYNLD HSHANYFCN LSIFDPPPFK VTLTGGYLHI YESQLCCQLK FWLPIGCAAF
VVVCILGCIL ICWLTKKKYS SSVHDPNGEY MFMRAVNTAK**

[REFERENCES]

- 1. Beier K.C., *et al.* (2000) *Eur. J. Immunol.* 30:3707-3717.**
- 2. Ling V., *et al.* (2001) *Genomics* 78:155-168.**
- 3. Aicher A., *et al.* (2000) *J. Immunol.* 164:4689-4696.**
- 4. Tezuka K., *et al.* (2000) *Biochem. Biophys. Res. Commun.* 276:335-345.**