



TEL:4006-871-227 Web:www.ybio.net Email:shybio@126.com

YBA139Hu01 100 μ g

Recombinant Tumor Necrosis Factor Related Apoptosis Inducing Ligand (TRAIL)

Organism Species: Homo sapiens (Human)

*Instruction
manual*

FOR IN VITRO USE AND RESEARCH USE ONLY
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

9th Edition (Revised in Jul, 2013)

[PROPERTIES]

Residues: Val114~Gly281

Tags: N-terminal His-Tag

Accession: Q6IBA9

Host: *E. coli*

Subcellular Location: Membrane; Single-pass

type II membrane protein.


Purity: >95%

Endotoxin Level: <1.0EU per 1 μ g (determined by the LAL method).

Formulation: Supplied as lyophilized form in PBS, pH7.4, containing 1mM DTT, 5% trehalose, 0.01% sarcosyl and preservative.

Predicted isoelectric point: 8.6

Predicted Molecular Mass: 21.0kDa

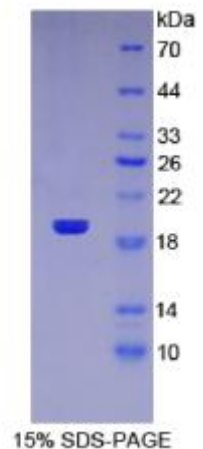



TEL:4006-871-227 Web:www.ybio.net Email:shybio@126.com

Applications: SDS-PAGE; WB; ELISA; IP.

(May be suitable for use in other assays to be determined by the end user.)

[USAGE]

Reconstitute in sterile PBS, pH7.2-pH7.4.



TEL:4006-871-227 Web:www.ybio.net

Email:shybio@126.com

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate of the target protein. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. (Referring from China Biological Products Standard, which was calculated by the Arrhenius equation.) The loss of this protein is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCES]

The sequence of the target protein is listed below.

VRERGPQ RVAAHITGTR GRSNTLSSPN SKNEKALGRK INSWESSRSG HSFLSNLHLR
NGELVIHEKG FYIYSQTYF RFQEEIKENT KNDKQMVQYI YKYTSYPDPI LLMKSARNSC
WSKDAEYGLY SIYQGGIFEL KENDRIFVSV TNEHLIDMDH EASFFGAFLV G

[REFERENCES]

1. Nguyen-Khac E., *et al.* (2010) Eur J Gastroenterol Hepatol 22:794-800.
2. Young R. P., *et al.* (2009) Postgrad Med J 85:515-524.
3. De Jager P. L., *et al.* (2009) Nat. Genet. 41:776-782.
4. Mavri A., *et al.* (2007) Thromb. Haemost. 97:62-66.