



TEL:4006-871-227 Web:www.ybio.net Email:shybio@126.com

YBA098Hu01 50µg
Recombinant Matrix Metalloproteinase 10 (MMP10)
Organism Species: Homo sapiens (Human)
Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

9th Edition (Revised in Jul, 2013)

[PROPERTIES]

Residues: Phe99~Cys476 (Accession # P09238), with N-terminal His-Tag.

Host: *E. coli*

Subcellular Location: Secreted, extracellular space, extracellular matrix.

Purity: >95%

Endotoxin Level: <1.0EU per 1µg (determined by the LAL method).


Formulation: Supplied as lyophilized form in PBS, pH7.4, containing 5% sucrose, 0.01% sarcosyl.

Predicted isoelectric point: 5.5

Predicted Molecular Mass: 44.5kDa

Applications: SDS-PAGE; WB; ELISA; IP.

(May be suitable for use in other assays to be determined by the end user.)

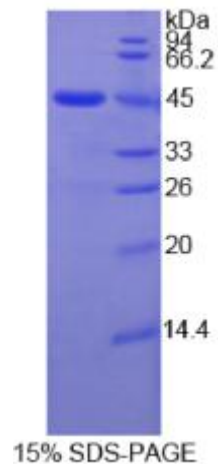

[USAGE]

Reconstitute in sterile PBS, pH7.2-pH7.4.



[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate of the target protein. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. (Referring from China Biological Products Standard, which was calculated by the Arrhenius equation.) The loss of this protein is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCES]

The target protein is fused with N-terminal His-Tag, its sequence is listed below.

MGHHHHHSGSEF-FS SFPMPKWRK THLYRIVNY TPDLPDAVD SAIEKALKVW
EEVTPLTFSR LYEGEADIMI SFAVKEHGDF YSFDGPGHSL AHAYPPGPGL YGDIHFDDDE
KWTEASGNTN LFLVAHELGL HSLGLFHSAN TEALMYPLYN SFTELAQFRL SQDDVNGIQS
LYGPPPASTE EPLVPTKSVP SGSEMPAKCD PALSFDAIST LRGEYLFFKD RYFWRRSHWN
PEPEFHLISA FWPSLPSYLD AAYEVNSRDT VFIFKGNEFW AIRGNEVQAG YPRGIHTLGF
PPTIRKIDAA VSDKEKKKTY FFAADKYWRF DENSQSMEQG FPRLIADDFP GVEPKVDAVL
QAFGFFYFFS GSSQFEFDPN ARMVTHILKS NSWLHC

[REFERENCES]

1. Bertini I., *et al.* (2004) *J. Mol. Biol.* 336:707-716.
2. Wertz K., *et al.* (2004) *Free Radic. Biol. Med.* 37:654-670.
3. Schnabl B., *et al.* (2005) *Biochem. Biophys. Res. Commun.* 335:87-96.
4. Aung P.P., *et al.* (2006) *Oncogene* 25:2546-2557.